


# belulu® ~classy~

## User's manual

Thank you very much for purchasing this product.

Please read this user's manual carefully before using this product.

There is a guarantee certificate at the end of this document,  
so please make sure to store it in a safe place.

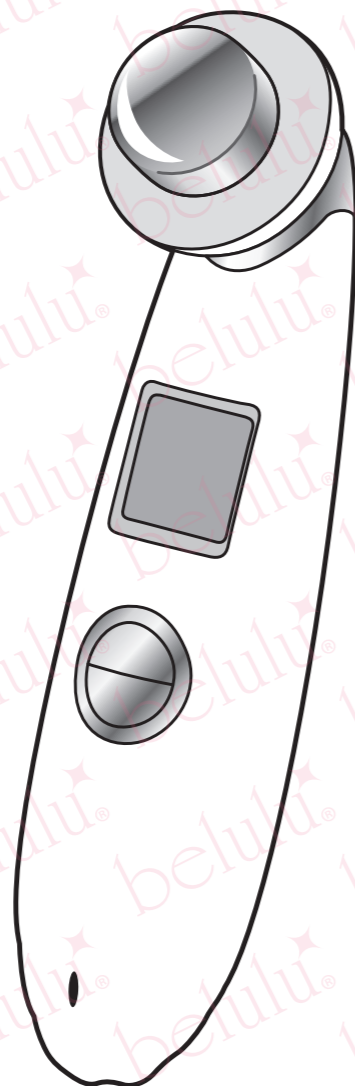

Please check the notice.

### ⚠ Notice

- Make sure to use the correct voltage, AC adapter, etc. by following the specifications in this document.
- When using this product, do not look directly into the LED light, or point it toward your eyes. Looking directly into the LED may be harmful to the eyes.
- Do not use this product if you find any problems with it, including any damages, faults, deformities, broken cables, etc.
- If you notice any smoke, offensive smell, abnormal heat, etc. coming from the product, turn off the power and stop using it immediately.
- Do not pull on the power adapter cable or touch the plug area with wet hands.
- Do not leave or use this product and the power adapter in places that may be subject to high temperatures, such as stoves, etc.
- Avoid using this product with wet hands, or when it is damp or covered in dirt, etc.
- This product is not waterproof. Never use this product in places with high humidity, such as bathrooms, etc. This may cause burns, electric shocks, or product faults.
- If you notice any abnormal effect on your skin, stop using this product and consult your physician.
- Look at the list below carefully, and do not use this product if any of those conditions apply to you.
- Never attempt to disassemble or reassemble or do any type of maintenance work on this product by yourself. This may cause the product or its accessories to become damaged or faulty. Any issues caused to the product or its accessories by these actions are not covered by the product guarantee, and the company shall bear no responsibility for them.

### ⚠ Caution

- Avoid using this product excessively on sensitive areas of your skin.
- If you feel any pain, etc. during use, move the 3D head away from your skin immediately.
- Do not subject the product to any strong actions such as by dropping or hitting it.
- The product is fully charged within approximately 4 to 8 hours. Do not continue to charge it once it is fully charged.
- Disconnect the charger cable from the product when it is not being charged.
- Do not wash this product with water. If dirty, take a damp cloth, wring it well, and use it to clean the product gently. Do not use solvents such as benzene, thinners, etc. to clean this product.
- Any injuries, faults, damages, etc. caused to the product due to your failure to follow the important notes and warnings stated above, any incorrect use or disassembly/reassembly of the product, or any maintenance work done by yourself, are not covered by the product guarantee, and the company shall bear no responsibility for them.

### 〈Do not use people who applies to the following〉


- Acute diseases ● Febrile illnesses ● Blood pressure disorders ● Infectious diseases
- Heart diseases ● Usage of medical electronic equipment (implanted medical devices, such as pacemakers, etc.; electronic life-support equipment, such as heart-lung machines, etc.; or wearable medical electronic devices, such as wearable electrocardiogram devices, etc.)
- History of skin cancer ● Skin disorders caused by diabetes or other diseases ● Currently pregnant ● Currently breast-feeding ● Skin conditions (atopic dermatitis, skin inflammation due to overexposure to the sun, infected and inflamed acne, and sensitive skin due to allergic reactions)
- Dermatitis caused by cosmetics ● Generally not in good physical health ● Usage of long-term doses of steroid hormones drugs, etc.

### 〈Do not use for these parts.〉


- Malignant tumors ● Cuts ● On bones or joints ● Areas with silicone due to plastic surgery, etc.
- Vocal cords area
- Delicate zones such as hair, around the eyes and ears, mucous membranes, etc.

## 2 Set / Name


※Accessories and packaging may be have changes

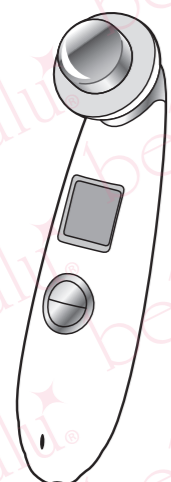

belulu Classy  
main unit

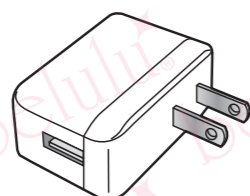


AC Adapter

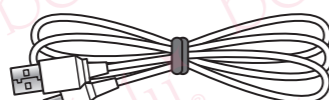


Charge cable



Light cut glasses  
(It is optional.)

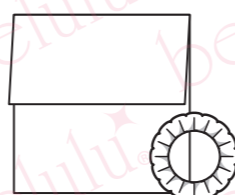

Fixing ring  
and facial cotton

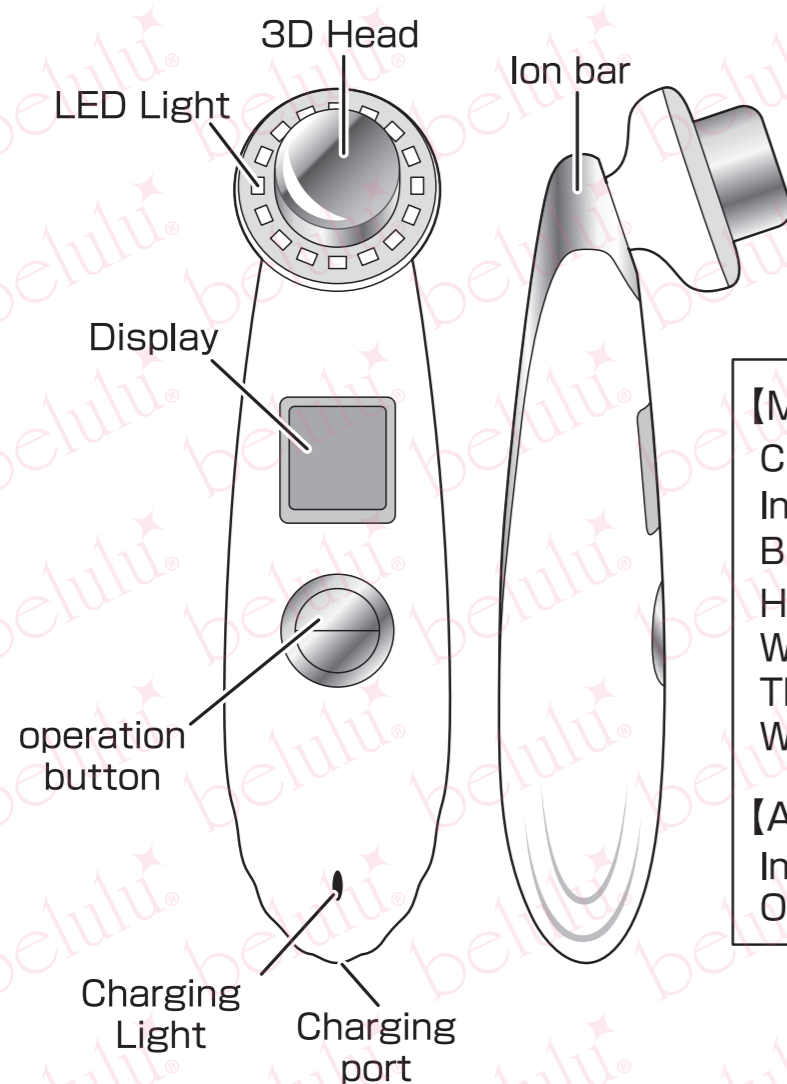

### 【Main Unit Specification】

Charging: 4~8 hours  
Input: 5V/0.5A  
Battery: Lithium-ion-battery  
Height: 160mm  
Width: 40mm(head)  
Thickness: 32mm(Body)  
Weight: 142g

### 【Adapter Specification】

Input: 100~240V 50/60Hz 0.2A  
Output: 5V/1A

## 3 Charging

●The charge notification is lower right in the display. Please charge when the remaining battery is 1 or 0.



full charge



Available



Available  
(Charging required)



Charging  
required

※Depending on the production period, the remaining battery level may be displayed immediately after a full charge, but this is due to a built-in program and there is no problem in use.

●Make sure to recharge the product in the following situations: if you are using the product for the first time after purchase, if you have not used it for a very long time, if its power seems low, or if it seems to work for short periods of time.

●Please turn the power off when you charge.

●Insert the charging cable's connector into the device's charging port and insert the AC adapter into a household power outlet.

●The charge notification light blinks with a red light.

※The device may become slightly hotter while charging; however, this is normal.

●The charging time and operating time will vary depending on how frequently the product is used.


●The charging light goes off when the product is fully charged. Make sure to disconnect the cable once charging is complete.

●This product cannot be used while connected to the charging cable. Make sure to check that the charging cable is disconnected when using the product.

■The most suitable room temperature for charging is 15°C to 30°C. Charging the product at low temperatures may shorten its battery life, so make sure to charge the product at the appropriate temperature. In addition, the product's charging speed may decrease at approximately 0°C or less.


■In order to protect the built-in battery, make sure to charge the product once per month even when you are not using it for long periods of time.

■If the product's operating time seems short even though it has been charged correctly, the cause of this is likely to be the battery life.

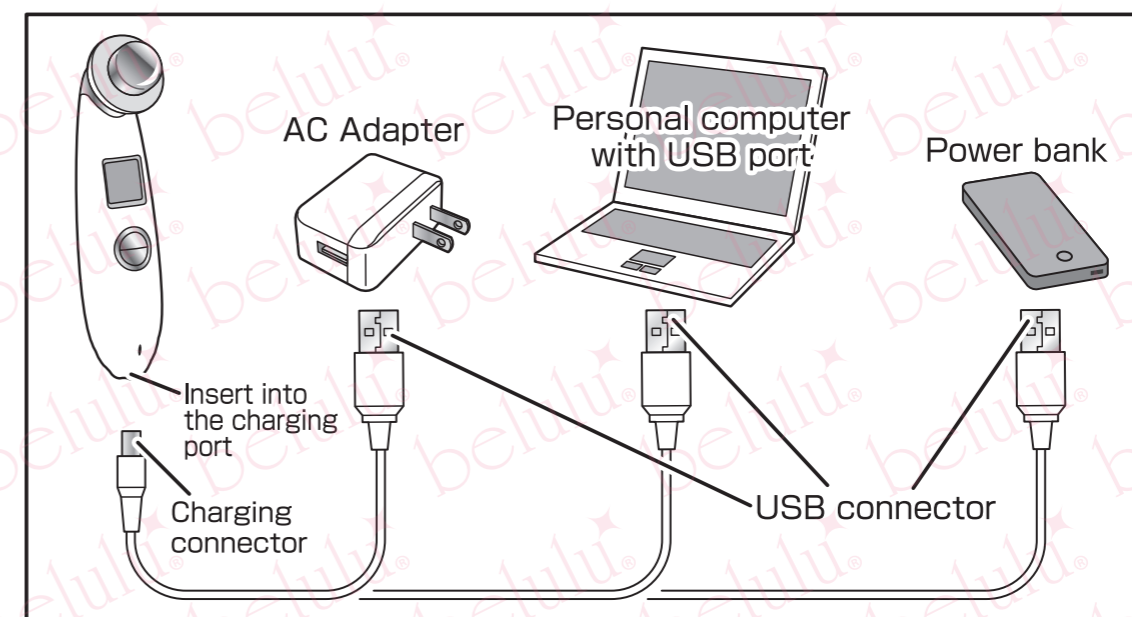

## 4 Operation

### ON/OFF button

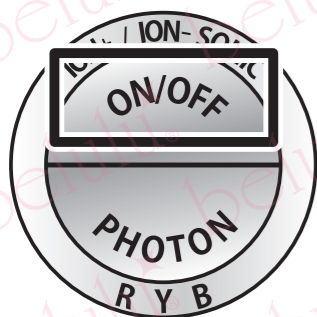


Please press the ON/OFF button to turn on the ion and sonic mode.


- Press the ON/OFF button over two seconds, the display will turn on and standby state.
- Then press the ON/OFF button each time the mode will change, Ion leading-out → Ion introduction → Vibration massage → Ion introduction + Sonic low → Ion introduction + Sonic medium → Ion introduction + Sonic high

After a certain time, the display light will be turned off. You cannot change the mode by pressing the button one time, it will just the display light on. When you change the mode, press the button one more time, after the display light turned on.

### Power on (Standby state)

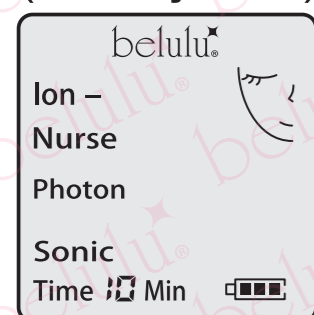


### Ion leading-out

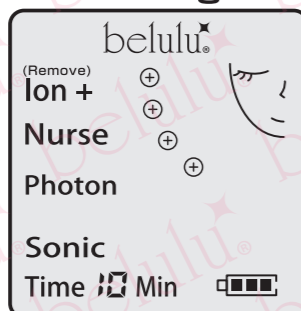


### Ion introduction Vibration massage




### Ion introduction Sonic [high]

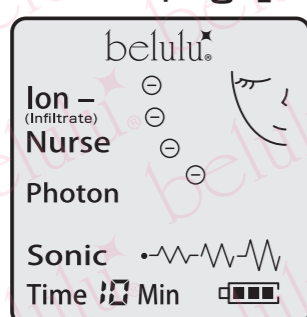

### Ion introduction Sonic [medium]

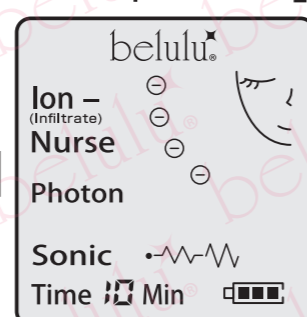


### Ion introduction Sonic [low]

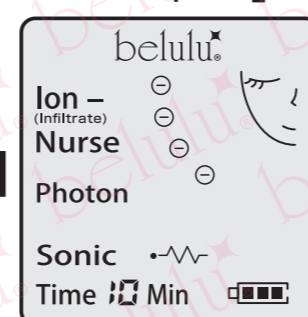

### PHOTON button


Choose the color of LED lights.

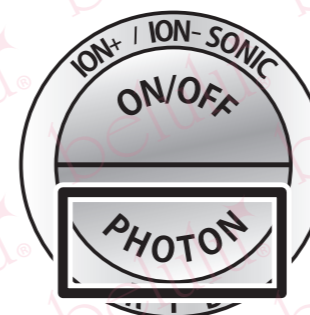


- Press the ON/OFF button over two seconds, the display will turn on and standby state.
- Then press the PHOTON button each time the light will change; Red LED → Yellow LED → Blue LED → LED turn off.
- ※ Please do not watch the LED light directly for a long time.

After a certain time, the display light will be turned off. You cannot change the LED light by pressing the button one time, it will just the display light on. When you change the LED light, press the button one more time, after the display light turned on.


### Power on (Standby state)

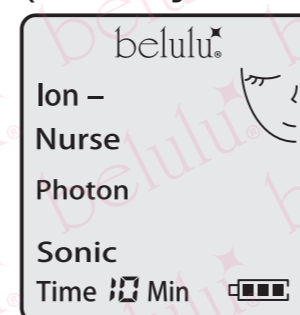


### Red LED

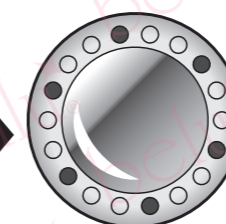

### Yellow LED

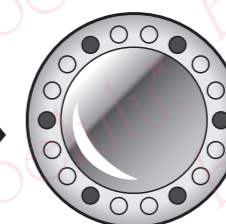

### Blue LED

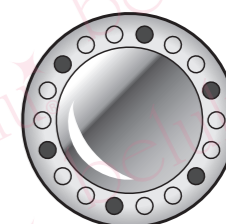

The mark of '☀' will be turn on during LED lightning.

You can use the Ion mode or Sonic mode at the same time by pressing the ON/OFF button.  
For more details of Ion mode or Sonic mode, please refer to 'ON/OFF button'.

You can use the LED photon at the same time by pressing the PHOTON button.

For more details of LED photo, please refer to 'PHOTON button'.

## 5 Before you use

Make sure to do an operation check and patch test before you use.

### Operation check

#### 【Ion leading-out】

- 1.Please refer to the 'ON/OFF button' and set the Ion leading-out mode.
- 2.You can see the '⊕' mark and it's normally.  
※There is no sounds or vibration during the Ion leading-out working.




#### 【Ion introduction】

- 1.Please refer to the 'ON/OFF button' and set the Ion introduction mode.
- 2.You can see the '⊖' mark and it's normally.
- 3.You can use the vibration massage and sonic at the same time.

#### 【Sonic】

- 1.Please refer to the 'ON/OFF button' and check the sonic mode switch; low → medium → high. The Ion introduction works at the same time.

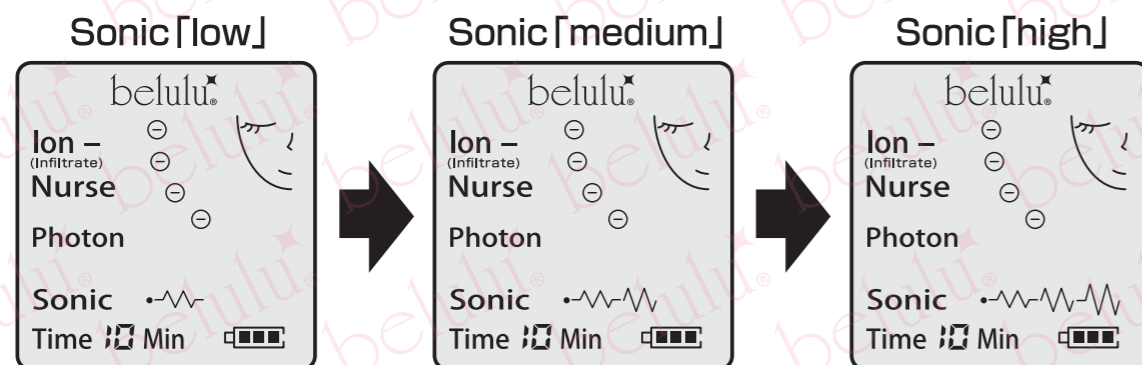

- 2.You can hear the small buzz sounds(it sounds like buzz). Please check the all, low, medium and high.
- 3.Put the small amount of the water on the 3D head, you can see the water waving or making some mist.  
※Please put a few drops of water.

#### 【LED photo】

- 1.Please refer to the 'PHOTON button' and check the LED light switch; Red → Yellow → Blue.
- 2.Please check the all the bulb working.  
※Please do not watch the LED light directly for a long time.
- 3.You can see the ☀ mark and it's normally.

### Patch test

- Wash and wipe your inside of arm.
- Apply some facial lotion to a weak part of your skin, such as the inside of your arm, and try out all the functionality provided by the product as outlined on the next page.
- Once you have verified that you have had no adverse reactions, you can move onto your face care. If you do experience any abnormality on your skin, do not use the product, and make sure to consult our company or your physician.

## 6 Function and effect

### 【Sonic】

High frequency stimulates the skin and improve the blood flow, lymph and metabolism.

Moisturizing, tightening up and good for improving wrinkle and dullness. It raises an effect of Ion introduction.

Ion introduction works at the same time and you cannot use separately.

### 【Ion leading-out】

For expelling pore dirt, grease off to deep clean skin, make the skin to rejuvenate.

### 【Ion introduction】

Can make the nutrition and moist penetrate into the deep skin for effective absorption and moist replenishment, accelerate skincare absorption

※Vibration massage works at the same time. You can't use separately.

### 【Vibration massage】

Vibration can make the nutrition and moist penetrate to the deep skin for effective absorption and moist replenishment, accelerate skincare absorption. Also help to improve sluggish blood circulation.

※Ion introduction works at the same time. You can't use separately.

### 【LED photon】

The care by LED light.

Each LED light color has different benefits.

You can change the color as you like.

#### Red LED

630nm

Will stimulate the production of Collagen is an essential protein used to repair damaged tissue and to replace old tissue. Best for fine lines and will also help to reduce large pores.

#### Blue LED

415nm(±10nm)

Bacteria in pimples contain Porphyrins; able to kill Porphyrins due to the high affinity. The blue light furthermore has calming action which is very effective for hypersensitivity.

#### Yellow LED

590nm(±10nm)

Can balance color pigment, reduce fine lines, nutrient aging skin, speed up the healing process of the wound, and lighten the scare.

## Combination

You can use some modes at the same time for time saving.

Ion leading-out

+

LED photon

Ion introduction(Vibration massage)

+

LED photon

Ion introduction(Sonic)

+

LED photon

## Count down timer

Power on and the 10 minutes timer stars. It turns off after 10 minutes automatically.


Please use the timer as a guide.

Select the mode and it indicates 'Time 10 Min'. It counts down every minute.



## 7 Let's start!

### Basic care

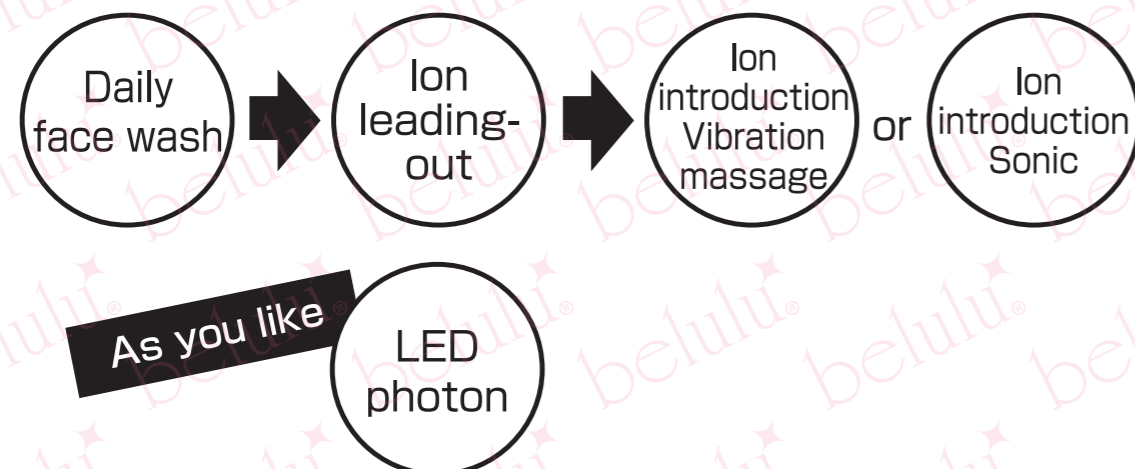

### Caution

- This product is not waterproof. Never use this product in places with high humidity, such as bathrooms, etc.
- Do not use this product in places that may be subject to high temperatures, such as stoves, etc.
- Do not watch the LED light directly for a long time.
- When doing your own skin care, make sure to use a gel or beauty cream. If you were to use this product on dry skin, it could cause the 3D head to move less smoothly, potentially damaging your skin. No need to apply the beauty serum, when you use only LED photon care.
- Grip the Ion bar.
- Apply the head area onto the skin by going down in a straight line. Also, do not press it down or rub it hard along the skin. Make sure to move the head gently and smoothly along the skin.
- Remove the all watch and jewelries
- If you notice any abnormality on your skin while using the product, or if you feel that it does not suit your skin in any way, stop using it immediately, and contact us or consult your physician.


### What to prepare

- Water or cleansing lotion  
...For Ion leading-out
- Beauty serum, gel or facial lotion  
...For Ion introduction. We recommend thick beauty serum which contains active ingredients blended, especially 'belulu b2 bottle'.
- Milky lotion or beauty cream  
...For all cares

**1. Ion leading-out** For 10 minutes / day, Three times / week

**Remove your makeup and wipe off moisture.**

Remove the dirt that you can't clean by ordinary face wash.


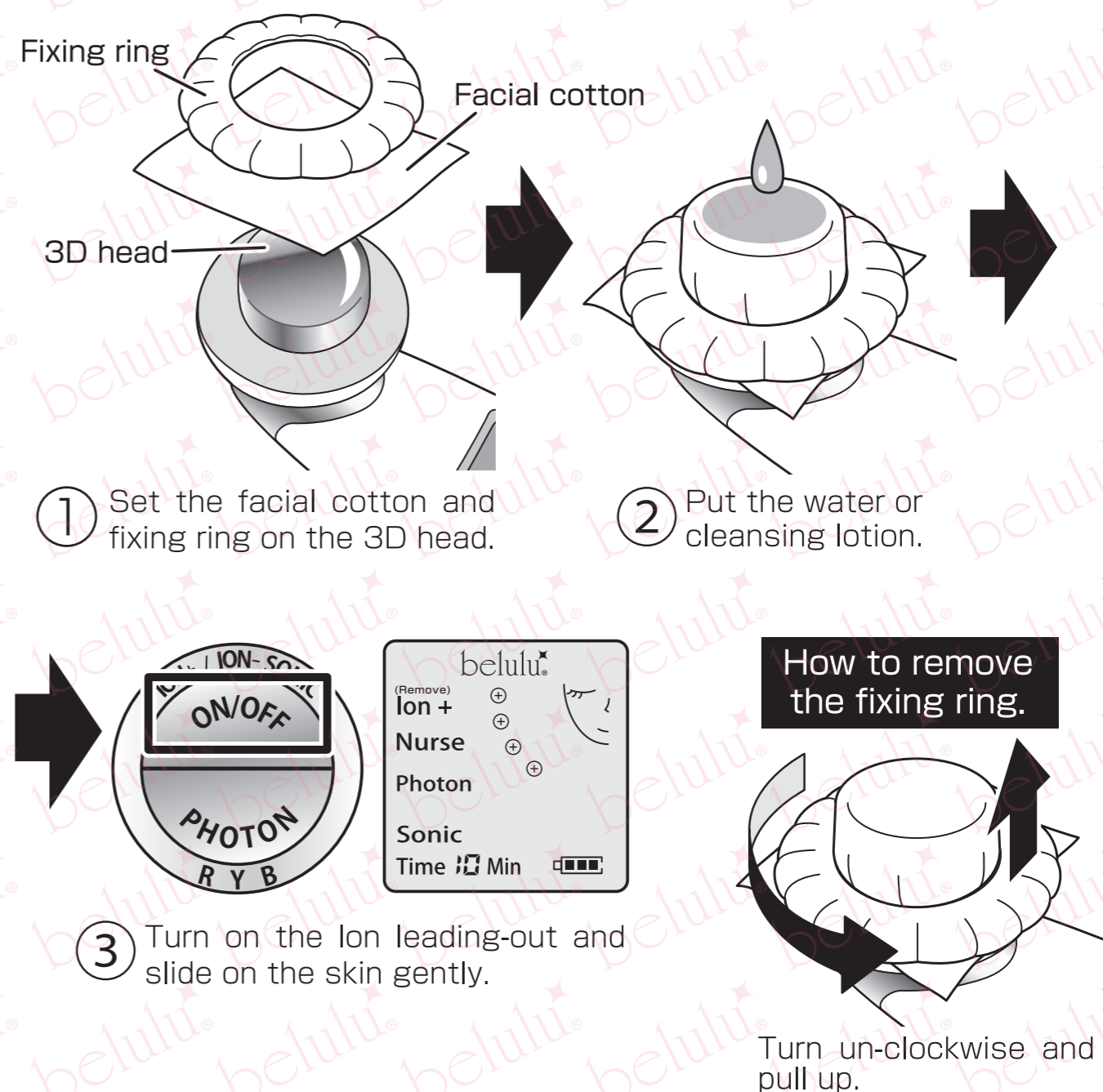


### Q&A Can I apply a facial lotion when using Ion leading-out?

Ion leading out is cleansing function. If you apply a facial lotion the pores will close.  
Apply water or cleansing lotion.

## [Point]


Slide the 3D head to each part as below image.  
Take a long care where you worry the dirt, especially around the nose.  
We recommend to add water or cleaning lotion for each time.

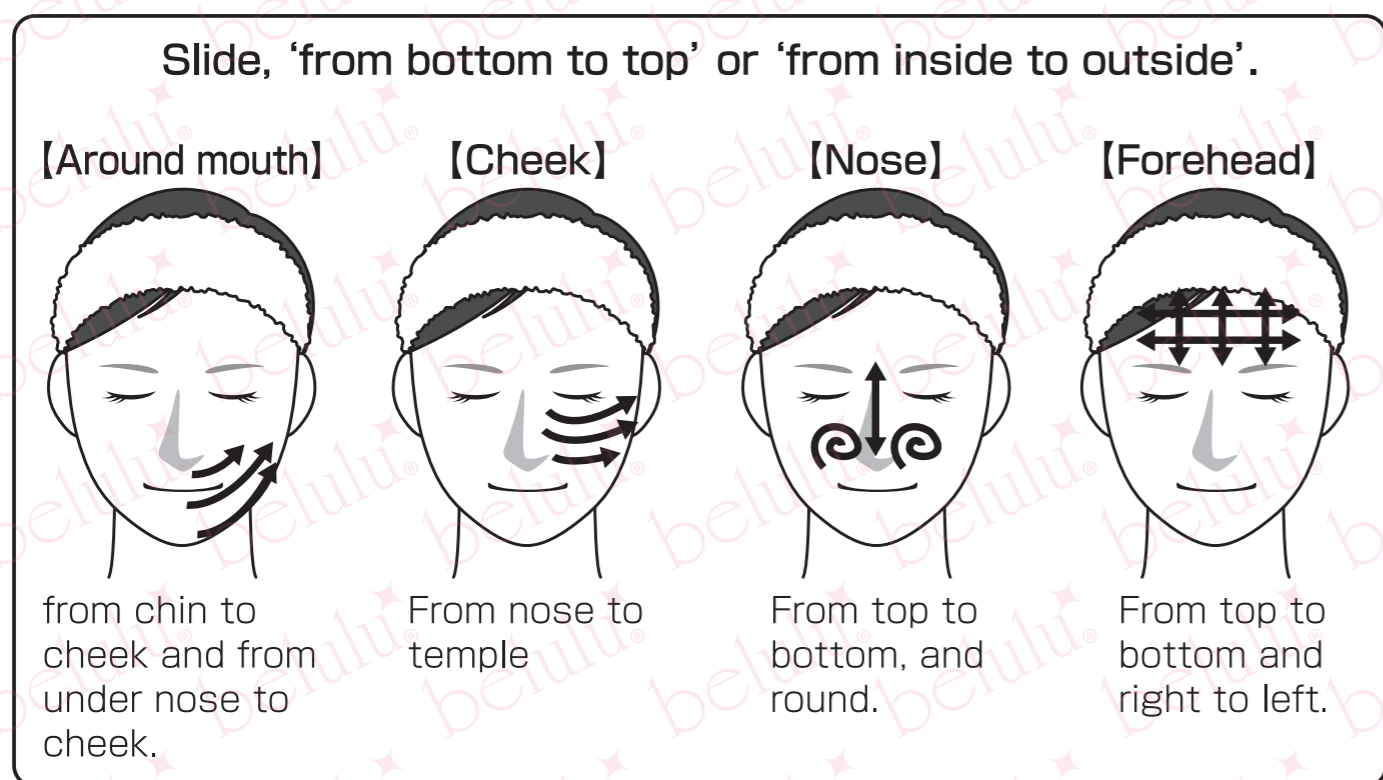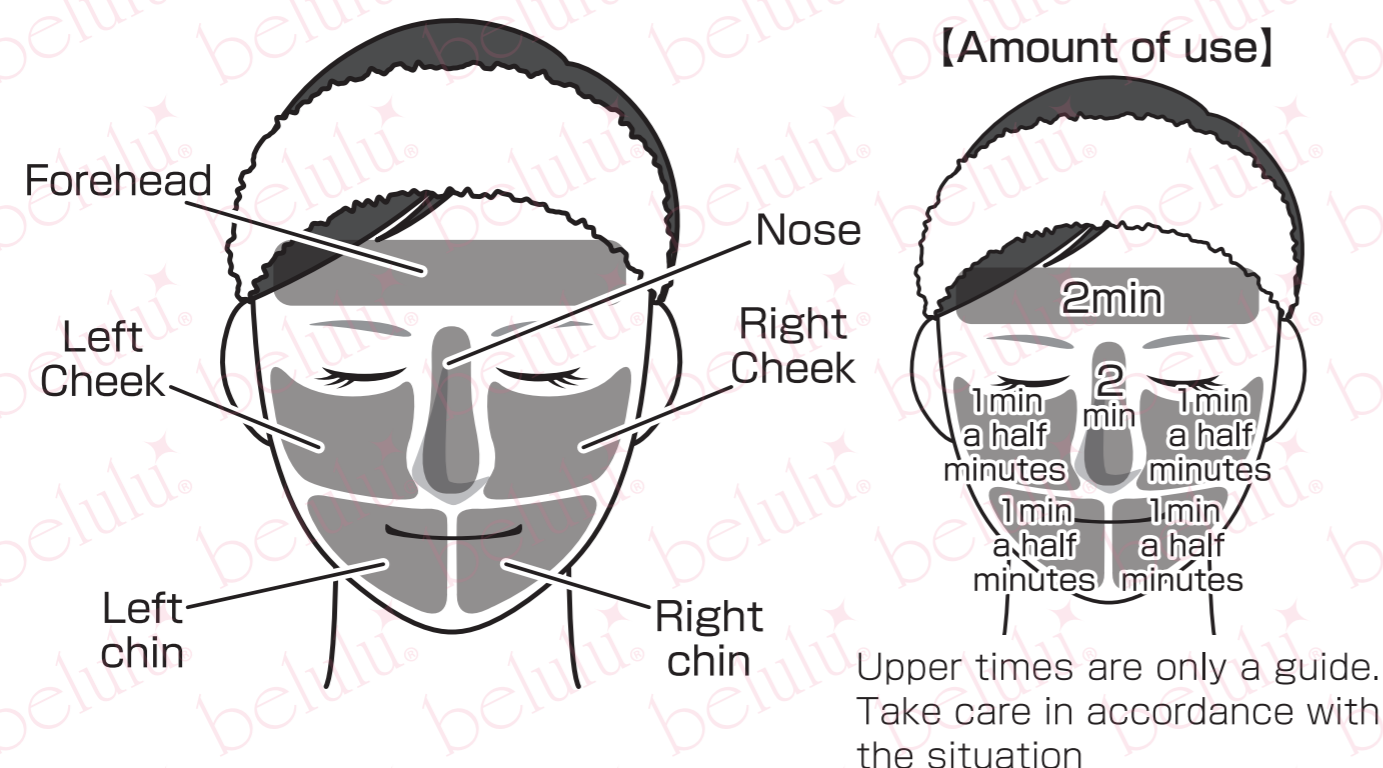



## 2-A.Ion introduction + Vibration massage 10 minutes / day, Everyday

### Wipe off moisture after Ion leading-out care.

Although it would not permeate with normal patting, it deeply moisturizes with negative ions produced from a weak electric current, allowing its beautifying elements to soak through. Choose this if you want to reduce swelling and have a brighter complexion from the vibrating massage that promotes good blood circulation.

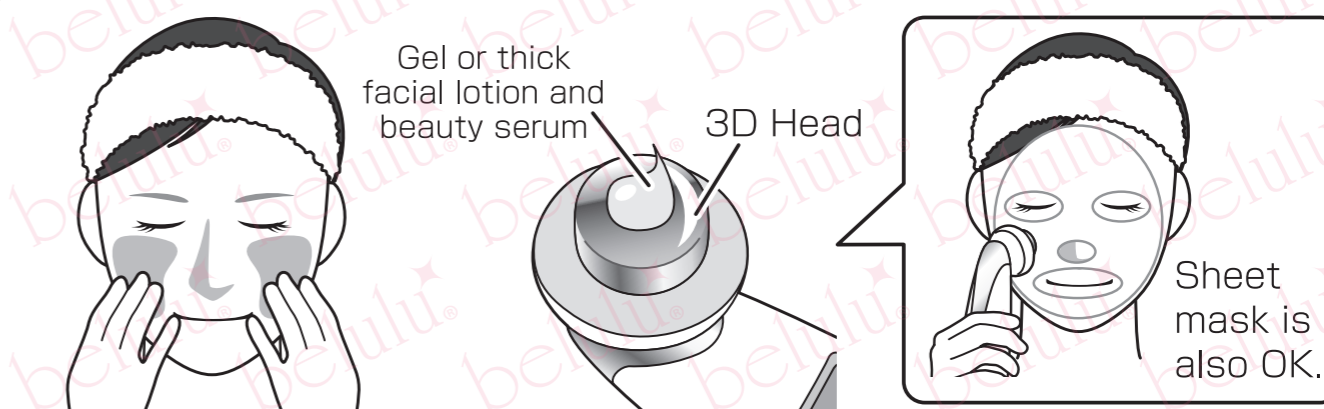


① Apply gel or thick facial lotion, beauty serum on your face. Or put on the 3D head and slide on your face.※Don't apply milky lotion or beauty cream.

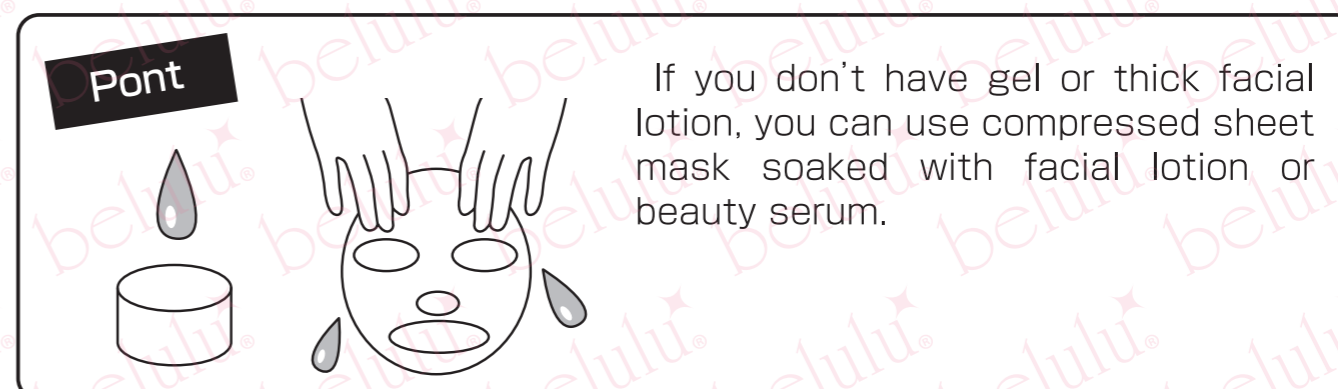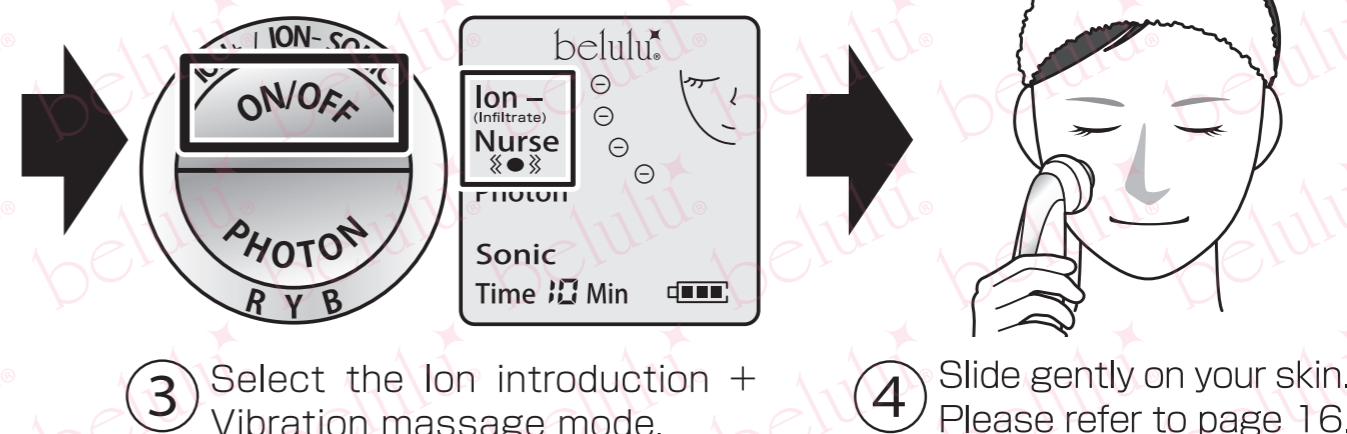


## 2-B.Ion introduction + Sonic 10 minutes / day, Three times / week

**Wipe off moisture after Ion leading-out care.**


Although it would not permeate with normal patting, it deeply moisturizes with negative ions produced from a weak electric current, allowing its beautifying elements to soak through. Choose this if you want to reduce wrinkles and have a dry skin from the sonic that promotes good blood circulation.

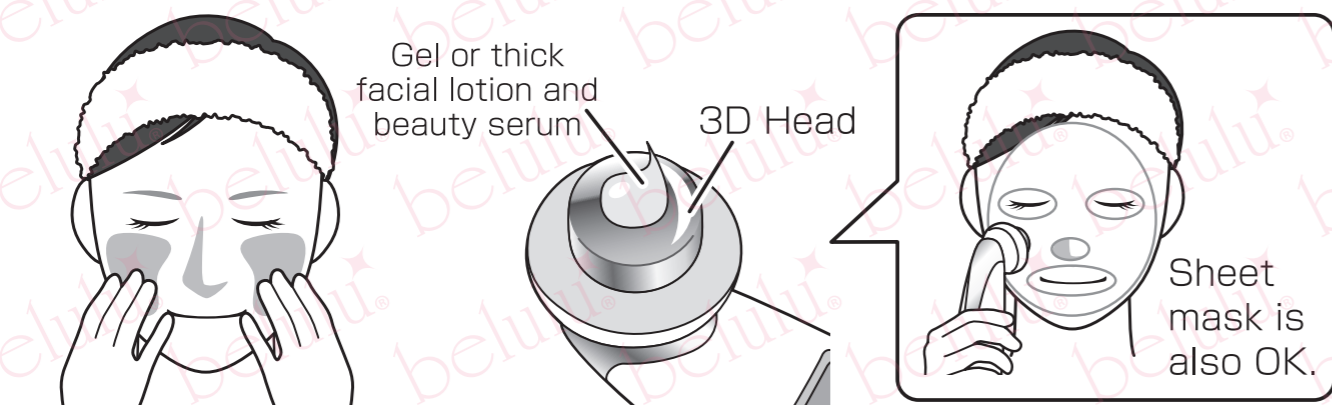

- ① Apply gel or thick facial lotion, beauty serum on your face. Or put on the 3D head and slide on your face.※Don't apply milky lotion or beauty cream.

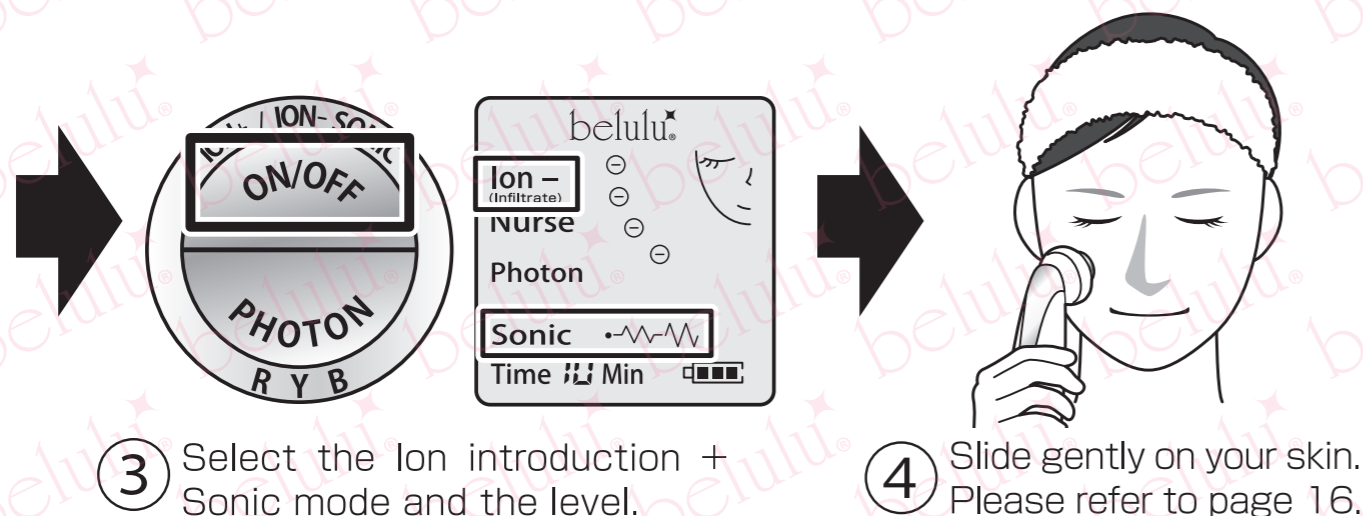


- ③ Select the Ion introduction + Sonic mode and the level.
- ④ Slide gently on your skin. Please refer to page 16.

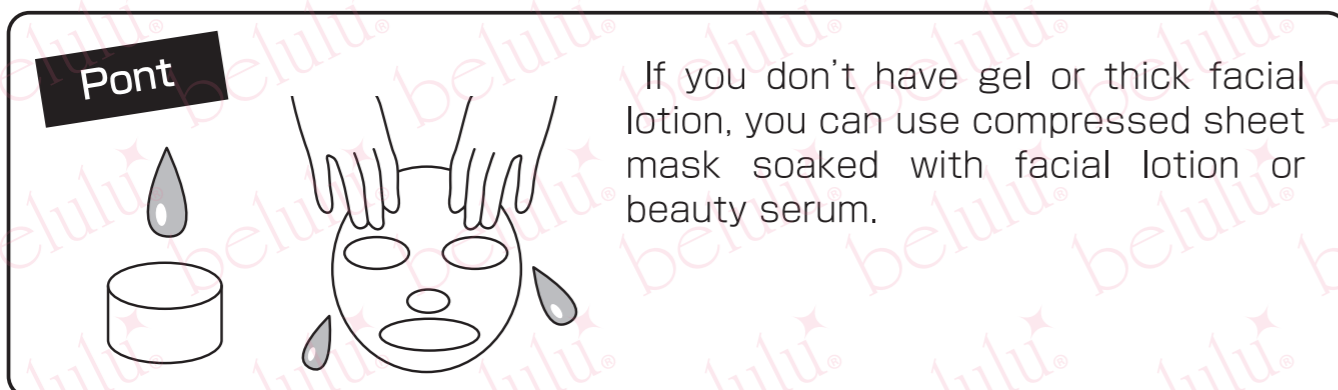

## 【Point】

Slide the 3D head to each part as below image.  
Take a long care where you worry dryness, especially around the mouth.  
We recommend to add gel or facial lotion, beauty serum for each time.

## 【お手入れ時間の目安】

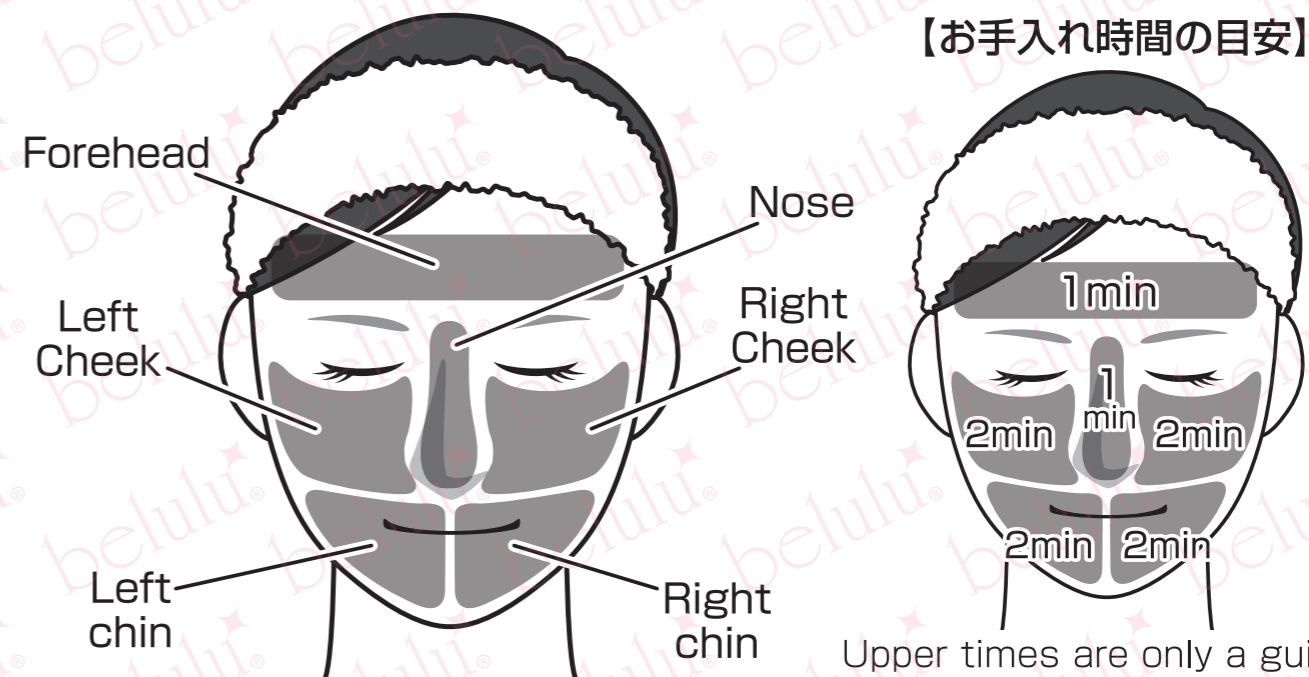

Upper times are only a guide.  
Take care in accordance with the situation

Slide, 'from bottom to top' or 'from inside to outside'.  
Slide the same as Ion leading-out care.

## 【Around mouth】



from chin to cheek and from under nose to cheek.

## 【Cheek】



From nose to temple

## 【Nose】



From top to bottom, and round.

## 【Forehead】



From top to bottom and right to left.

If this is the point where you are finishing up with your care routine, carefully blend in the rest of the gel, facial lotion, or beauty serum on your face, and then seal it with a milky lotion or beauty cream to compensate for the fat content.

LED photon Each color 10 minute/ day Everyday


After Ion introduction, leave and use the rest of gel, facial lotion, or beauty serum on your face.

Remove the sheet mask.

※Don't apply milky lotion or beauty cream.

In case of only LED photon care, you don't need apply gel.

Do not watch the LED light directly for a long time and wear light cut glasses.

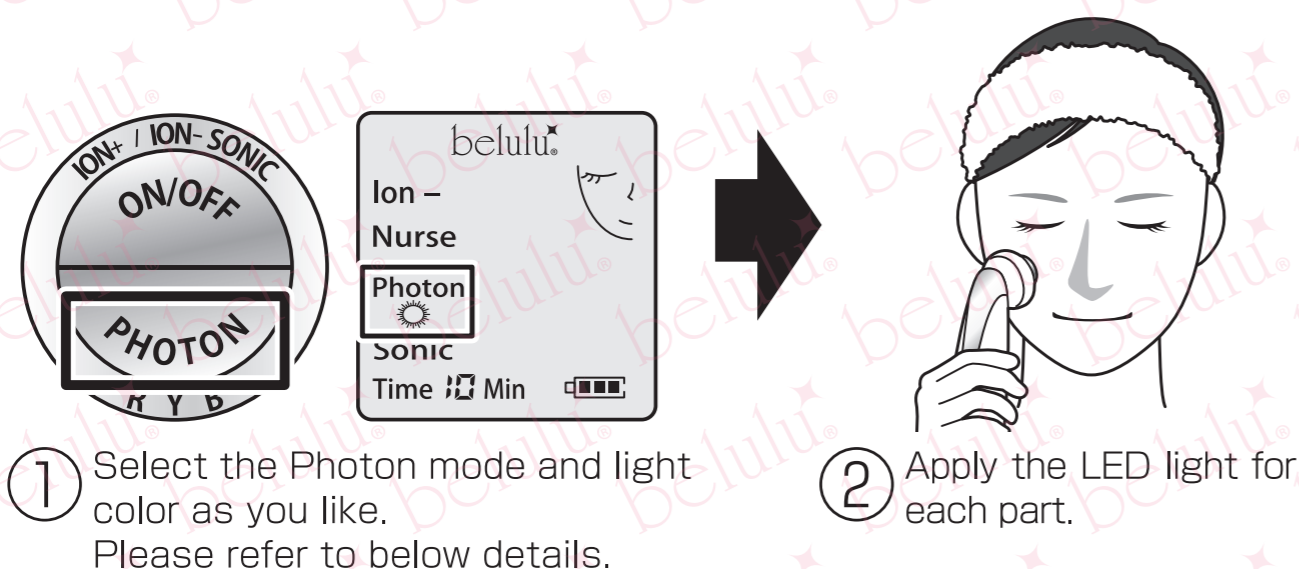

**Red LED**

630nm

Will stimulate the production of Collagen is an essential protein used to repair damaged tissue and to replace old tissue. Best for fine lines and will also help to reduce large pores.

**Blue LED**

415nm(±10nm)

Bacteria in pimples contain Porphyrins; able to kill Porphyrins due to the high affinity. The blue light furthermore has calming action which is very effective for hypersensitivity.

**Yellow LED**

590nm(±10nm)

Can balance color pigment, reduce fine lines, nutrient aging skin, speed up the healing process of the wound, and lighten the scare.

## 8 After care

- Turn off the power.
- Blend in the rest of the gel, facial lotion, or beauty serum on your face.
- Seal it with a milky lotion or beauty cream to compensate for the fat content.
- Wipe the gel or facial lotion, beauty serum on the 3D head and keep it clean.
  - ※ Do not wash this product with water. If dirty, take a damp cloth, wring it well, and use it to clean the product gently.
  - ※ Do not use solvents such as benzine, thinners, etc. to clean this product.
- Keep it in a box and avoid dust.
- Avoid storing in places with high temperatures, high humidity, or under direct sunlight, as well as places that are highly dusty or within the reach of children.

9

Q&A

The device doesn't work.	<p>●<b>It may be out of charge.</b> The rechargeable battery has the peculiar attribute of losing power a little at a time even when not turned on, and there may be times when it has no charge. Recharge for one hour and try again.</p> <p>●<b>Press the ON/OFF button over two seconds, the display will turn on and standby state.</b></p> <p><b>【Still have problem?】</b> We have taken early failures and breakdowns of this product into consideration, so please contact us.</p>
Not charged, not stable charging.	<p>●<b>The cable is connected securely?</b> It can't be charged if it isn't energized. Please refer the 'Charging' and check the charging port and connect securely.</p> <p><b>【Still have problem?】</b> We have taken early failures and breakdowns of this product, cable or battery into consideration, so please contact us.</p>
The sonic doesn't work.	<p>●<b>Check the operation.</b> The ultrasonic waves produce roughly 3 million oscillations per second, and the oscillations at that speed are neither visible nor are they felt by the skin. Please refer to '5. Before you use' and 'Sonic' and check the operation.</p> <p><b>【Still have problem?】</b> We have taken early failures and breakdowns of this product into consideration, so please contact us.</p>

Ion leading-out and Ion introduction don't work.	<p>●<b>Check the operation.</b> Due to Ion leading-out and Ion introduction's feature of sending a weak electric current into the skin, the oscillations or stimulus cannot be heard. ※You can hear the Vibration massage sounds. Please refer to '5. Before you use' and 'Ion leading-out', 'Ion introduction' and check the operation.</p> <p><b>【Still have problem?】</b> We have taken early failures and breakdowns of this product into consideration, so please contact us.</p>
The display light turned off	After a certain time, the display light will be turned off. Push the any button one more time.
The mode doesn't change.	After a certain time, the display light will be turned off. You cannot change the mode by pressing the button one time, it will just the display light on. When you change the mode, press the button one more time, after the display light turned on.
Display is not visible and LED light is weak.	<p>●<b>Low charge</b> If its power seems low, please recharge.</p> <p><b>【Still have problem?】</b> We have taken early failures and breakdowns of this product into consideration, so please contact us.</p>
During use, there is a "beep" and it powers off.	It powers off after 10 minutes automatically.

**When in use, my hands will tingle as if it is short circuiting.**

●**This is not a short circuit.**

A tiny amount of electricity is flowing through the metal handle of this device, and it is meant to pass through your skin's barrier to permeate its beautifying elements deep in the skin, so it is possible depending on the person that this tiny electric current will be felt. This small electric current is harmless to your body, but if you are concerned, use it by turning down the strength level of the ultrasonic waves.

**【Still have problem?】**

Please contact us. Within 7 days from arrival.

## Q&A

**Q. Must I use the beauty serum?**

**A.** We recommend thick beauty serum which contains active ingredients blended, especially 'belulu b2 bottle'.

**Q. Is it safe to use if I'm pregnant?**

**A.** You can use it even if you are pregnant, but your skin is in a delicate condition from the hormonal imbalance during pregnancy, so depending on the timing when it is used, it is possible in rare cases that it will harm your skin. It is recommended that you make sure to consult your physician and use it while observing how it affects your skin.  
(The beauty device's functions or the ingredients from the beauty serum have no effect on a mother's body or a fetus, so please don't worry about that.)

**Q. What is the 3D head made from?**

**A.** It is stainless.  
It is strong against cracks, corrosion, or rusting, and it is hard for it to trigger allergies. It may not be difficult for it to trigger allergies for everyone.  
Please do a patch test.

**Q. What should I do if it isn't suitable for my skin?**

**A.** We will refund within 7 days from arrival, subject to shipping fee.

**Q. My skin is stiff, so right after I wash my face, I want to use Ion leading-out after applying the facial lotion.**

**A.** Ion leading-out has the function of absorbing the dirt not removed from washing your face, but if you apply facial lotion right after washing your face, the dirt will remain and clog your pores, so there will be no effect from the Ion leading-out and Ion introduction after that. Stiff skin is caused by chronic dryness of the inner part of the skin, so you will need to thoroughly moisturize it. It will hurt at first, but strictly adhere to carefully removing the dirt by Ion leading-out and then moisturize by Ion introduction. Continual use moisturizes the inner part of the skin, improving stiff or dry skin.

**Q. Even after using it, I can't feel any of its effects.**

**A.** How the effects materialize is different depending on the person. This device was manufactured with the intention of it being used at home, so it does not have the same power as from a beauty salon, and there is no immediate effect. By continuing to use it every day, you will start to see its effects.

**Q. Even after Ion leading-out, the cotton is not brown.**

**A.** How the dirt is removed is different depending on the person, so the cotton will not always necessarily turn brown. The brown stain is from foundation, but for anyone who does not normally wear that much makeup, or if the foundation is cleanly removed by cleansing and washing of the face, the cotton will not turn brown. Blocked pores and oils can also be removed by drawing out the ions, but as they are white, most of that will assimilate into the white cotton making it hard to see.

**Q. Should I take time every day to perform this care routine?**

**A.** This device's built in functions demonstrate highly beautifying effects, but overuse will put a burden on your skin. So as not to burden your skin during use, abide by the frequency of use as listed in this manual. Frequent use every day will not produce immediate effects, it is important to carefully continue to use it on a daily basis.

If you have any questions or concerns,  
please feel free to contact us.

**JAPAN: [support@kireido.com](mailto:support@kireido.com)**