

belulu Triple Rich Booster

User's manual

Thank you very much for purchasing this product. Please read this user's manual carefully before using this product. There is a guarantee certificate at the end of this document, so please make sure to store it in a safe place.



PLEASE READ BEFORE USING


Warning

- Make sure to use the correct voltage, AC adapter, etc. by following the specifications in this document.
- Do not use this product if you find any problems with it, including any damages, faults, deformities, broken cables, etc.
- If you notice any smoke, offensive smell, abnormal heat, etc. coming from the product, turn off the power and stop using it immediately.
- Do not pull on the power adapter cable or touch the plug area with wet hands.
- Do not leave or use this product and the power adapter in places that may be subject to high temperatures, such as stoves, etc.
- Avoid using this product with wet hands, or when it is damp or covered in dirt, etc.
- This product is not waterproof. Never use this product in places with high humidity, such as bathrooms, etc. This may cause burns, electric shocks, or product faults.
- If you notice any abnormal effect on your skin, stop using this product and consult your physician.
- **Never attempt to disassemble or reassemble or do any type of maintenance work on this product by yourself. This may cause the product or its accessories to become damaged or faulty. Any issues caused to the product or its accessories by these actions are not covered by the product guarantee, and the company shall bear no responsibility for them.**

Caution

- Do not too much use this product for the same area.
- If you notice any abnormal effect on your skin, stop using this product.
- Do not subject the product to any strong actions such as by dropping or hitting it.
- Do not wash this product with water. If dirty, take a damp cloth, wring it well, and use it to clean the product gently. Do not use solvents such as benzine, thinners, etc. to clean this product.
- **Any injuries, faults, damages, etc. caused to the product due to your failure to follow the important notes and warnings stated above, any incorrect use or disassembly/reassembly of the product, or any maintenance work done by yourself, are not covered by the product guarantee, and the company shall bear no responsibility for them.**

Spec

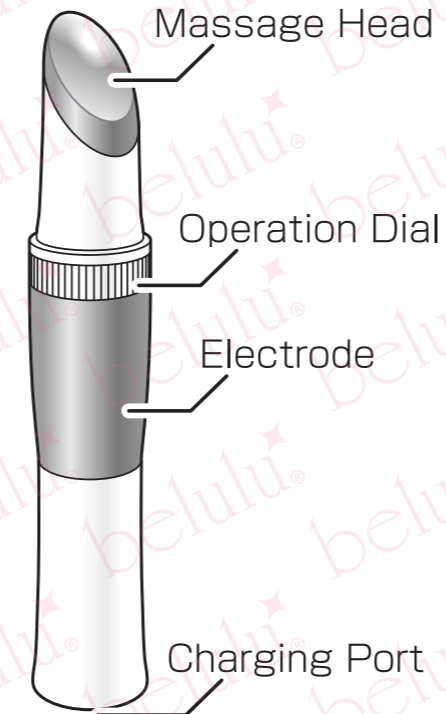


Charging Time:
2-6 Hours

Size:
Diameter: 20 mm
Height: 128 mm
Weight: 160 g

【Specified USB Type
AC Adapter】



Input: 100~240V
50/60Hz 0.2A
Output: 5V/1A


How to use

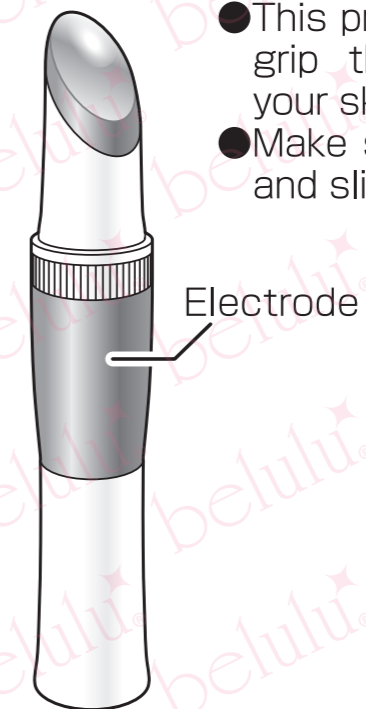

- This product operates when you grip the electrode and touch your skin at the same time.
- Make sure to grip the electrode and slide on your skin.

- Make sure to recharge the product in the following situations: if you are using the product for the first time after purchase, if you have not used it for a very long time, if its power seems low, or if it seems to work for short periods of time.
- Insert the charging cable's connector into the device's charging port and insert the AC adapter into a household power outlet.
- The charge notification light blinks with a blue light.
- ※ The device may become slightly hotter while charging; however, this is normal.
- The charging time and operating time will vary depending on how frequently the product is used.
- The charging light goes off when the product is fully charged. Make sure to disconnect the cable once charging is complete.
- This product cannot be used while connected to the charging cable. Make sure to check that the charging cable is disconnected when using the product.
- The most suitable room temperature for charging is 15°C to 30°C. Charging the product at low temperatures may shorten its battery life, so make sure to charge the product at the appropriate temperature. In addition, the product's charging speed may decrease at approximately 0°C or less.
- In order to protect the built-in battery, make sure to charge the product once per month even when you are not using it for long periods of time.
- If the product's operating time seems short even though it has been charged correctly, the cause of this is likely to be the battery life.

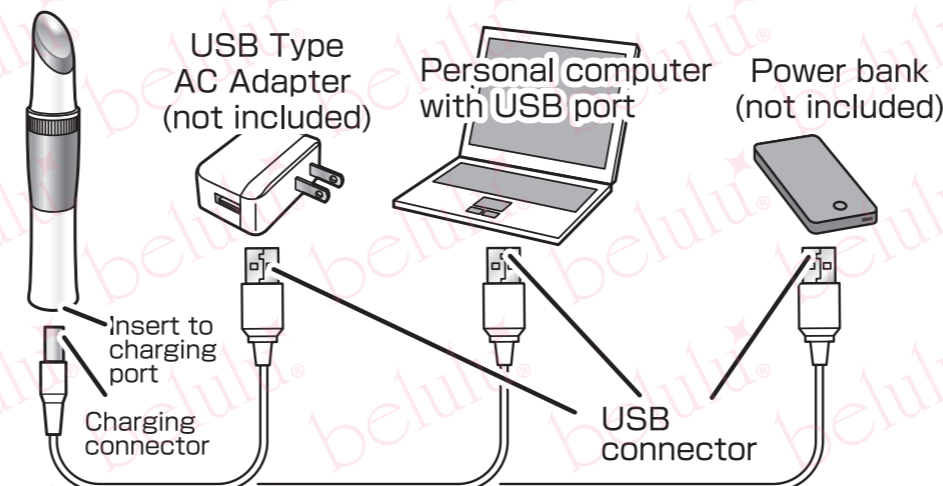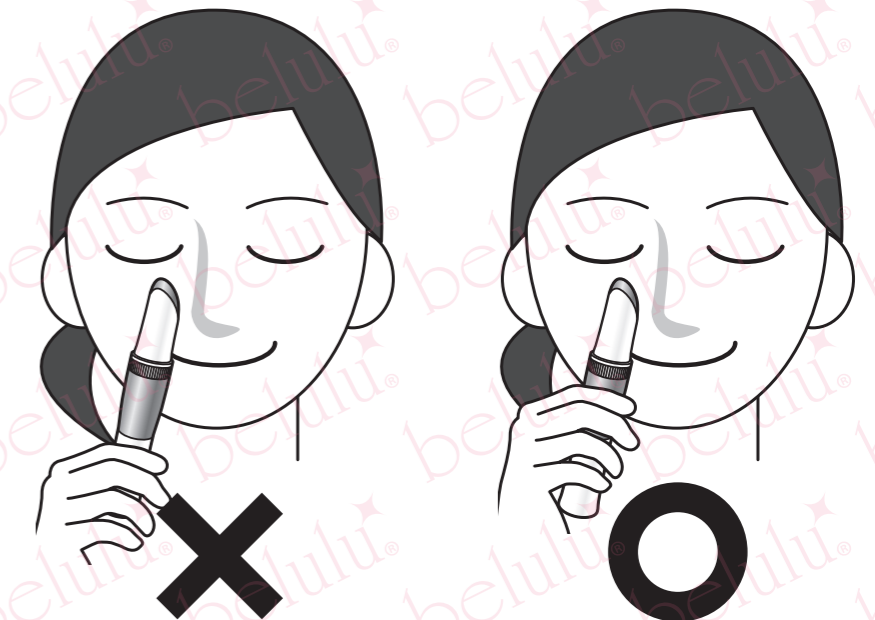


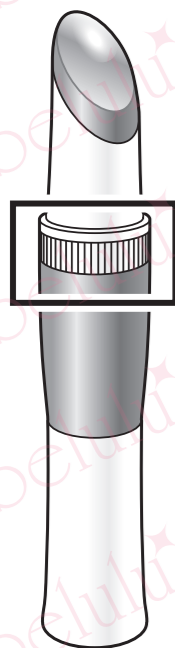

You can switch Vibration mode or Vibration + Heating mode by setting the Operation Dial.

After care

- Blend in the rest of the gel or beauty serum on your face.
- It will stop by touching away from your skin.
- Wipe the gel on the head and keep it clean.
- ※ Do not wash this product with water.
- ※ Do not use solvents such as benzene, thinners, etc. to clean this product.
- Avoid storing in places with high temperatures, high humidity, or under direct sunlight, as well as places that are highly dusty or within the reach of children.

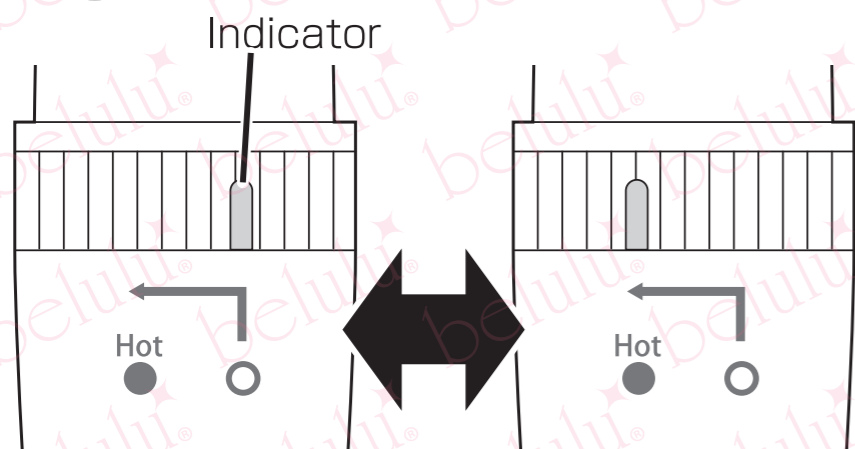

When you set the indicator on the mark of 'O', it will operate 'Vibration mode'.

When you set the indicator on the mark of 'Hot', it will operate 'Vibration + Heating mode'.



Please apply any beauty gel or beauty serum and slide gently around your eyes.

Do not apply on your eyeballs.

If you have any questions or concerns, please feel free to contact us.

JAPAN: info@kireido.com
CHINA: china@kireido.com